

Cincinnati Sub-Zero

Z-PLUS BENCHTOP

Temperature / Humidity Chambers



Engineered &
Manufactured
in the USA

www.cszproducts.com

Stackable Benchtop Model Saves Floor Space Expandable for Future Applications

Flexible to meet your test space and requirements

Our new [8 cu. ft. benchtop test chamber](#) is ideal for the most common testing applications, available as temperature only unit or temperature and humidity chamber. This particular model is designed for multiple use cases, offering flexibility and the ability to test more product in the same workspace. The chamber may be used in a variety of configurations:

1. Standard benchtop model, placed on existing table top.
2. Single chamber with optional cart provides storage for test equipment on the supplied shelf, saving space.
3. Stacked double chamber configuration with base saves floors space, **doubling test capacity in the same footprint. Run two different tests in separate chambers without taking up additional space in your lab.**



2 Single Unit with Optional Cart



Stacked, Double Unit with Base



Cabinet Features

The contemporary and smooth cabinet design is aesthetically pleasing to fit the style of your laboratory. Stainless steel interior with oven-baked, powder-coated exterior for added durability.

- Two, 4" access ports centered on the left-hand and right hand sides for ease of cable routing. Ports are fully welded to eliminate leaks and increase chamber life
- One adjustable product shelf slides out to provide easier access to your product. The shelf design is non tipping and supports large product loads
- Single-handed, lockable latch operation for ease of use
- Fog-free viewing window provides product viewing
- Interior light illuminates chamber workspace and product
- Removable side panels allow for easy access to all systems
- Double gasketing assures a vapor tight seal and insulation minimizes heat loss
- Each chamber is fully pressure tested to eliminate potential leaks



Z-Plus Benchtop Temperature/Humidity Chambers

Refrigeration Features

- Environmentally safe refrigerants are non-flammable, and have a zero Ozone Depletion Potential (ODP)
- Refrigeration service taps and refrigeration pressure gauges are included for easy maintenance
- Refrigeration system saves energy costs and prevents coil frost up for efficient operation


Our patented Tundra® system offers increased performance, reduces maintenance and energy usage, saving up to 54% on operational costs. Available in temperature ranges as low as -45°C. See pg. 7 for details.

Electrical & Safety Features

- All wiring between the electrical panel and the various components pass through terminal strips for easy maintenance. Wiring is color coded and identified with numbers that correspond to system schematics
- Multiple safety devices include a primary over-temperature limit to protect your product and an over-temperature safety limit to protect the chamber. All branch power and control circuits are individually fused
- UL-508A listed electrical panel

Communications

- Z-Plus models come equipped with TCP/IP and EIA-232 serial communications, for monitoring & control, and optional IEEE-488 GPIB
- Safety relay connection is provided to protect your device under test by removing power to it when the chamber is not running

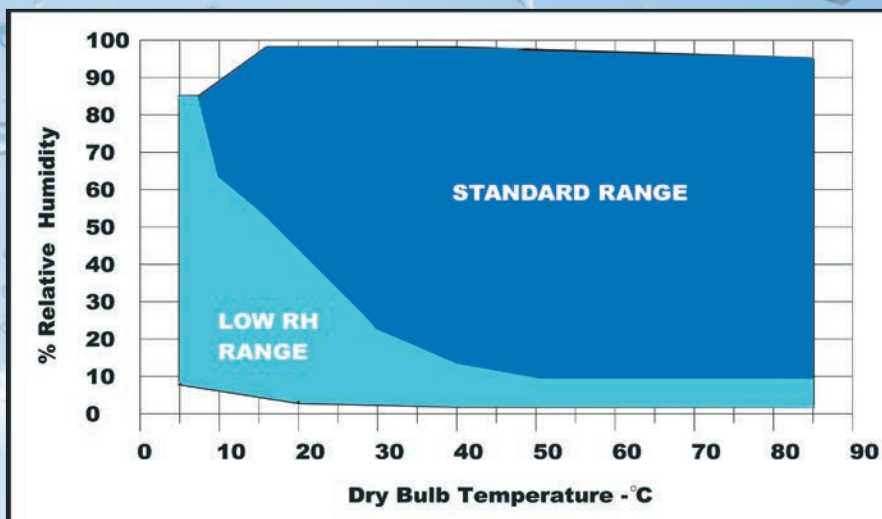

Humidification System (ZPBH Models)

- The Z-Plus features our fast response, tight control, humidity system. Humidification systems allow for humidity and moisture resistance testing meeting a variety of commercial and military standards
- These systems include an electronic, solid-state humidity sensor for accuracy and minimal maintenance



PC Software Controls Multiple CSZ Chambers

Using our optional EZ-View windows software, from one central PC you are able to monitor, control, datalog, create profiles, and receive alarm notifications. This software package allows you to control up to 31 chambers.



Choose from a complete selection of options for flexibility and variety of configurations.

Optional Accessories

- Dry Air Purge
- Recirculating Water Supply
- Humidity Water Filtration
- Demineralizer Cartridge Filters
- Low RH
- Water Pressure Regulator
- Running Time Meter
- IEEE-488 Computer Interface
- Windows-Based Software
- Chart Recorder
- Additional Shelves
- Customer Event Digital Outputs
- Digital Inputs
- Heat Only
- Blank Door without window
- Main Power Cord
- Main Power Disconnect
- Redundant Product Hi/Low Limit
- Reinforced Chamber Floor
- Reinforced Shelf
- EZ-View PC Software
- Single Unit Cart
- Stacked, Double Unit Base

Z-Plus Benchtop Performance Specifications

ZPB(H) - 8	
Workspace Volume	8 Cu Ft (230 L)
Workspace Dimensions WxDxH in (mm)	24"W x 24"D x 24"H (610 x 610 x 610)
Exterior Dimensions WxDxH in (mm)	36.5" x 80.5" x 36" (927 x 2044 x 914)
Exterior w/ Single Unit Cart* WxDxH in (mm)	36.5"W x 80.5"D x 74"H (927 x 2044 x 1879)
Exterior w/ Stacked & Base* WxDxH in (mm)	36.5"W x 80.5"D x 78"H (927 x 2044 x 1981)
Temperature Range	Tundra®: -45°C to +190°C (-49°F to 375°F)
Temperature Control Tolerance**	±0.5°C (±0.5°F) at temperatures >0°C ±1.1°C (±2°F) at temperatures <0°C
Humidity Range Optional Range	10% to 98% RH 5% to 98% RH
Humidity Control Tolerance	±5% RH at steady state conditions after stabilization
Distributed Shelf Load Capacity	110 lbs.

*Optional

**At steady state conditions after stabilization.

Model	Cooling Performance with Empty Chamber in Minutes from:				Cooling Rate ² °C / min	Heating Performance with Empty Chamber in Minutes from:		
	24°C		85°C	24°C		-34°C		
	-18°C	-34°C	-40°C	-40°C		93°C	190°C	24°C
ZPB-8-2-SCT/AC	6	12	18	30	4.1	10	30	10
ZPB-8-3.5-SCT/AC	4	8	12	22	5.6	5.5	16	5.5

Model	Live Load Capacity Watts			Electrical Power Requirements Full Load Amps ¹	
	-18°C	-34°C	-40°C	208/240V, 1 Phase	208/240V, 3 Phase
ZPB-8-2-SCT/AC	1800	950	725	31	27
ZPB-8-3.5-SCT/AC	2110	1100	925	-	37

Performance is based on 240V, 60 Hz. operation and a 24°C ambient. Lower voltages will extend heating performance.

¹ Electrical requirements based on temperature only units. Amperage may increase on humidity models. See quotation for actual values.

² From +85°C to -40°C




Go Green with our Tundra® Cooling System

Tundra is a patented refrigeration system design that uses one compressor and is available in temperature ranges as low as -45°C. This system provides additional performance for fast transition rates, increased live load capacity and can save up to 54% on energy costs.



The Tundra systems offers these benefits:

- Increased Performance
- High Reliability
- Quiet Operation
- Energy Efficient
- Improved Service ability
- Reduced Utility and Maintenance Cost
- Increased Value

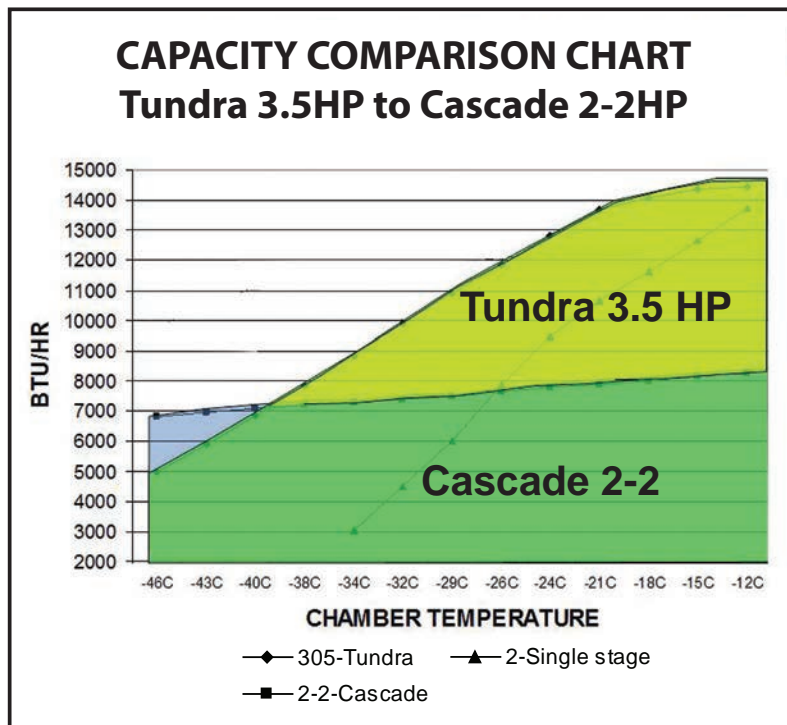

EZT-570 Touchscreen Controller Find out more

The Next Generation Controller with Smartphone Technology

All features are built into the controller interface so no additional software or internet is required for access to all the features the controller has to offer.

Communications & Connectivity

- Monitor and/or Control the chamber remotely for anytime, anywhere access from any device using LAN VNC
- Alarm notification sends email and/or text messages
- Email built-in to send data, alarm, audit trail files directly from controller
- Ethernet TCP/IP and serial communications



Save valuable time with the ease of use of the EZT-570 featuring fewer steps to accomplish your daily testing needs while incorporating simplified operation and programming to test faster.



Profiling

- Profiling includes up to 99 steps and 1000 cycles
- Program ramp steps entering time or °C/min
- Programs may be written using product control function
- Easily review profile using trend chart or review list of steps before running profile
- Profile status view displays current step, estimated start/stop date and time and more
- Profiles may be transferred to different chambers via USB or optional EZ-View software
- Automated delay profile start

Data Logging

- Configurable log interval, data file length, filename, operator entered batch & lot information as well as an unlimited number of operator notes saved to the data file
- Access data files directly from controller or PC
- Easily download profiles, alarm files, audit trail files and data files using USB or email from controller in a compatible .csv file format for ease of use. Also import profiles to other chambers saving valuable profile entry time
- Files may also be automatically backed up daily for hassle-free file management using FTP. FTP/FileWeb/DataWeb (LAN/WAN)





User Convenience & Flexibility

- Controller may be configured in 11 languages (additional available)
- Selectable power failure/recovery options
- Full system security allows up to 30 different users with four different levels of security
- Audit trail files track changes in settings by each user
- Configure alarm setting and maintenance alerts (optional)

Graphing Technology

- Real-time trend display graph with adjustable time and min/max values
- Four configurable trend graphs with left & right axis (additional available)
- Graph historical data files
- Zoom in/out of graphs for a closer look

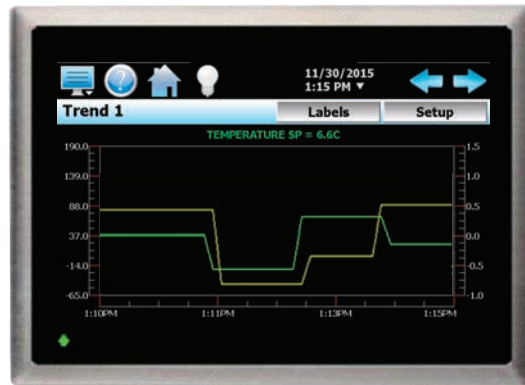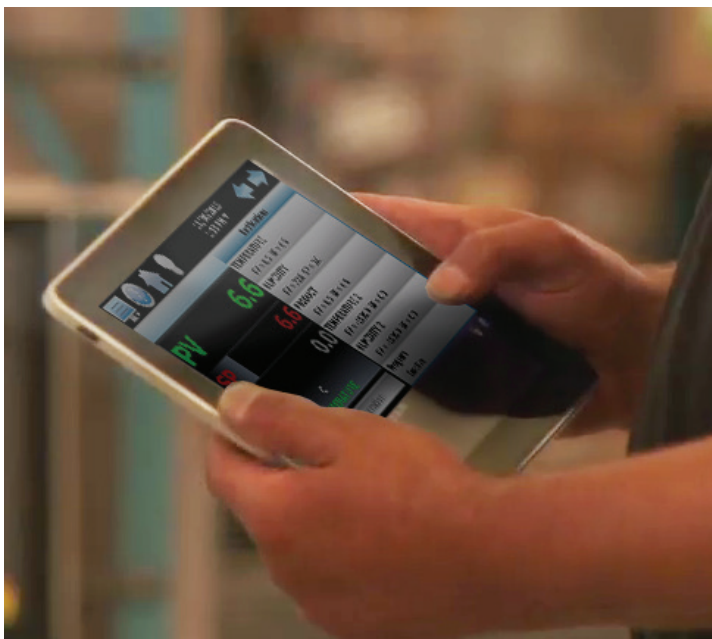

Enhanced Communications & Control Options

- Digital input option provides 8 inputs that can be configured for various control functions including starting, stopping and pausing a profile. "Wait for" function allows the user to pause a profile during a particular step of the profile until a specific digital input is turned on or off
- Digital output "customer event" feature provides 15 programmable outputs. Each output can be configured to perform other operations including alarm or profile status indicators for more control over your testing
- Optional refrigeration monitor package displays and data logs temperatures and refrigeration system compressor suction/discharge pressures
- Condensation control option helps prevent condensation from collecting on the part by automatically managing the air dewpoint
- Bar code option allows user to scan barcode to start profile and to add notes to current data file when datalogging



Cincinnati Sub-Zero is a product brand of Weiss Technik North America, Inc. Weiss Technik North America is a member of the Weiss Technik group of companies, a division of the Schunk Group with its headquarters in Heuchelheim, Germany. Weiss Technik is the world's largest manufacturer of environmental simulation systems and employs more than 2,900 people in 18 group companies in 15 countries.

[Request Information](#)



CSZ
Cincinnati Sub-Zero
The Testing Standard.

Weiss Technik North America, Inc.
Cincinnati Facility
12011 Mosteller Road
Cincinnati, OH 45241

(p) **513-326-5252**

(f) **513-326-5258**

www.cszproducts.com

